Supplying Hospitals  
for Over 90 Years

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## User Instructions

Safety Summary: The following recommendations are made for the safe assembly of the raised toilet seat.

### WARNINGS:


- ⚠ **DO NOT** assemble or use this equipment without first reading and understanding this instruction sheet. If you are unable to understand the Warnings, Cautions or Instructions, contact a healthcare professional, dealer or technical personnel before attempting to install this equipment. Otherwise, injury or damage may occur.
- ✓ Certain toilet fixtures will prevent the locking raised toilet seat from working properly and/or will damage the locking raised toilet seat. If the toilet bowl has on its top or inner ring any raised or protruding points or sections, the locking raised toilet seat should not be used.
- ✓ **BEFORE** each use, ensure that **ALL** attaching hardware is securely tightened.
- ✓ Use the arm supports **ONLY** for assistance. **DO NOT** attempt to use the arms to support full body weight.
- ✓ Always check to make sure that the raised toilet seat is properly and securely locked in place before using. Make certain that the knob that secures the raised toilet seat to the commode is tight at all times.
- ✓ Center your weight on the raised toilet seat. The seat may tip if you lean too far forward to either side, or if equipped by leaning on only one armrest. Evenly distribute weight between the armrests during use.
- ✓ The raised toilet seat's weight capacity is: 350 lbs. (159 kg.)
- ⚠ Users with limited physical capabilities should be supervised or assisted when using the raised toilet seat.
- ✓ Rotate the adjustment knob by hand only to prevent over-torque of the assembly. Use of a wrench or other tool will break the screw assembly. **DO NOT FORCE** knob as this may damage the knob and/or the raised toilet seat.
- ✓ The Medline raised toilet seat is not designed for permanent installation. The unit should be loosened and retightened at least once a week to ensure a proper fit.
- ⚠ If, for any reason, a proper fit in installation cannot be achieved, **DO NOT** use and contact retailer where purchased for return instructions.

### To Install the Equipment to the Toilet Bowl

The Medline raised toilet seat has been designed to accommodate a wide variety of toilet bowl types. However, certain bowl types will require the removal of the existing seat and/or lid to ensure proper fit when installed.

1. Lift the seat and lid of the toilet.
2. Rotate adjustment knob of the raised toilet seat counter clockwise until the paddle is far enough back to clear the inner rim of the bowl.
3. Place the raised toilet seat onto the bowl rim hooking the rear wings under edge of bowl. Check that the unit rests flat on the bowl rim. The knob must be at the front end of the toilet bowl.
4. Center the unit on the bowl rim and push the back as far as possible with your hands.
5. Rotate adjustment knob clockwise to tighten. Turn until resistance is felt, then turn an additional 1/8" until snug.

Rear Seat Wings
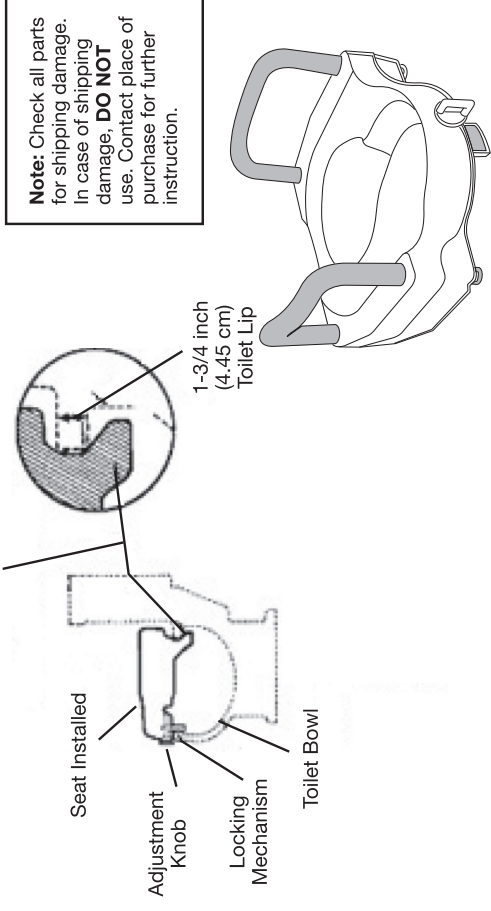
### To Use Toilet with this Equipment

1. Carefully back up to the raised toilet seat until you feel the back of your knees touching it.
2. Reach back for both armrests to guide you into position.
3. Slowly lower yourself onto the raised toilet seat while bending your elbows and knees.
4. Carefully use the armrests for assistance.
5. Slowly raise yourself off of the toilet seat.
6. Ensure that you have your balance before releasing your grasp from the raised toilet seat.

### To Remove the Equipment from the Toilet Bowl

1. Rotate the knob counterclockwise until seat is loose.
2. Lift up the front of the raised toilet seat until the adjustment bracket clears the lip of the toilet bowl.
3. To completely remove the seat from the toilet bowl, pull the seat towards you at an angle.

### To Clean the Equipment

This raised toilet seat can be easily cleaned with regular household cleaning products. "DO NOT" use abrasive powders or scouring pads.

1. Clean the seat at least once a week, using a non-abrasive detergent or non-abrasive cleaner with warm water.
2. Rinse with warm water and wipe with a non-abrasive cloth.

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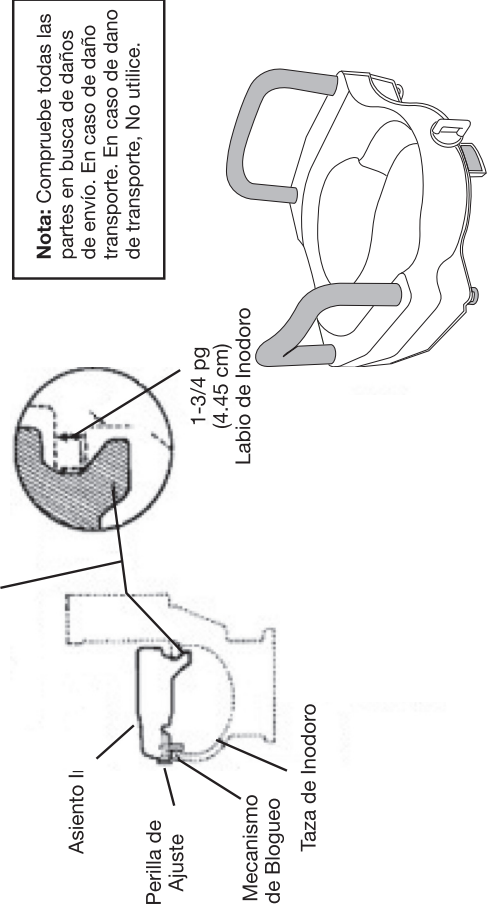
LIMITED WARRANTY

Your Medline product is warranted to be free of defects in materials and workmanship for the duration of use by the original purchaser. This device was built to precise standards and carefully inspected prior to shipment. This limited warranty is an expression of our confidence in the materials and workmanship of our products and assures years of convenient, dependable service. In the event of a defect not covered by this warranty, we will, at our option, repair or replace the device. This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, which are subject to normal wear and replacement.

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© 2014 Medline Industries, Inc.
Medline is a registered trademark of Medline Industries, Inc.
Medline Industries, Inc.
Medline es una marca registrada de Medline Industries, Inc.